
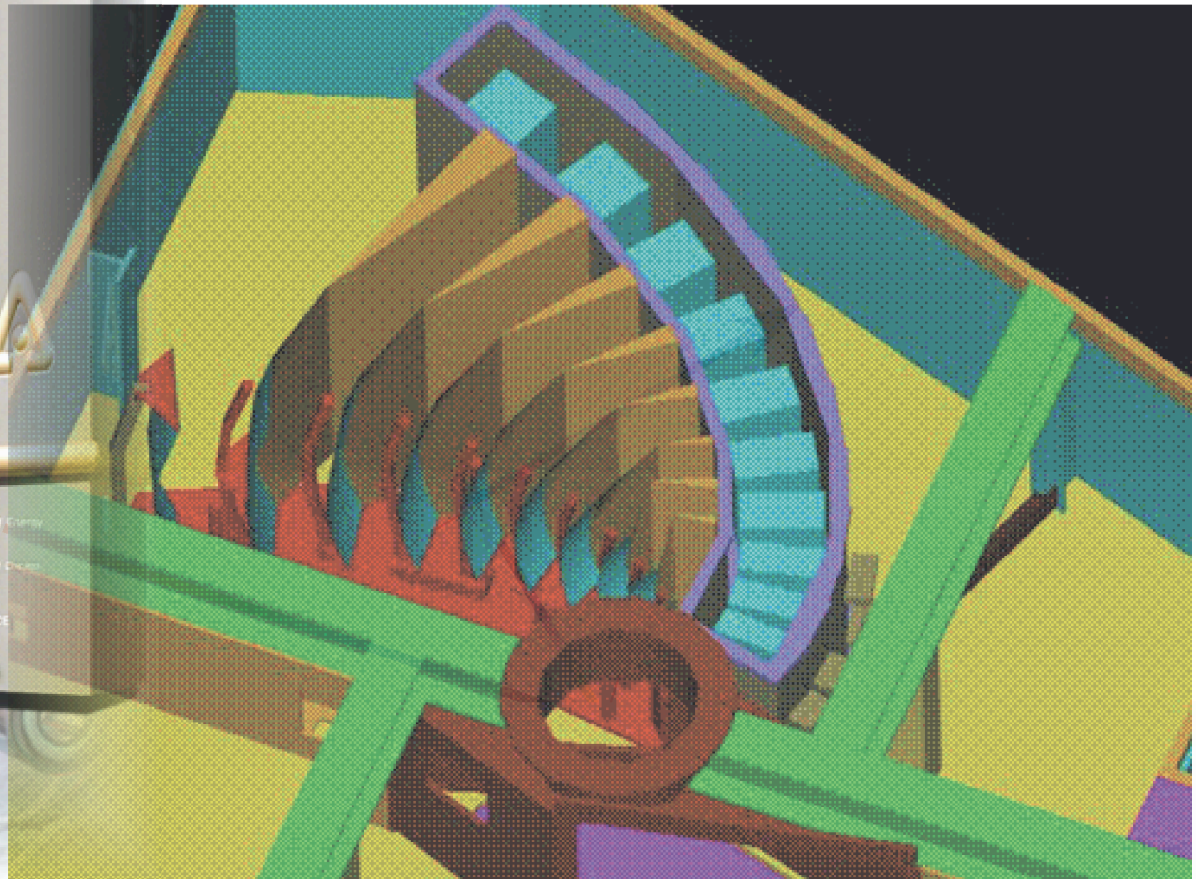


Intense Pulsed Neutron Source (IPNS)


<http://www.pns.anl.gov>

The QENS User Manual
Dec 2003

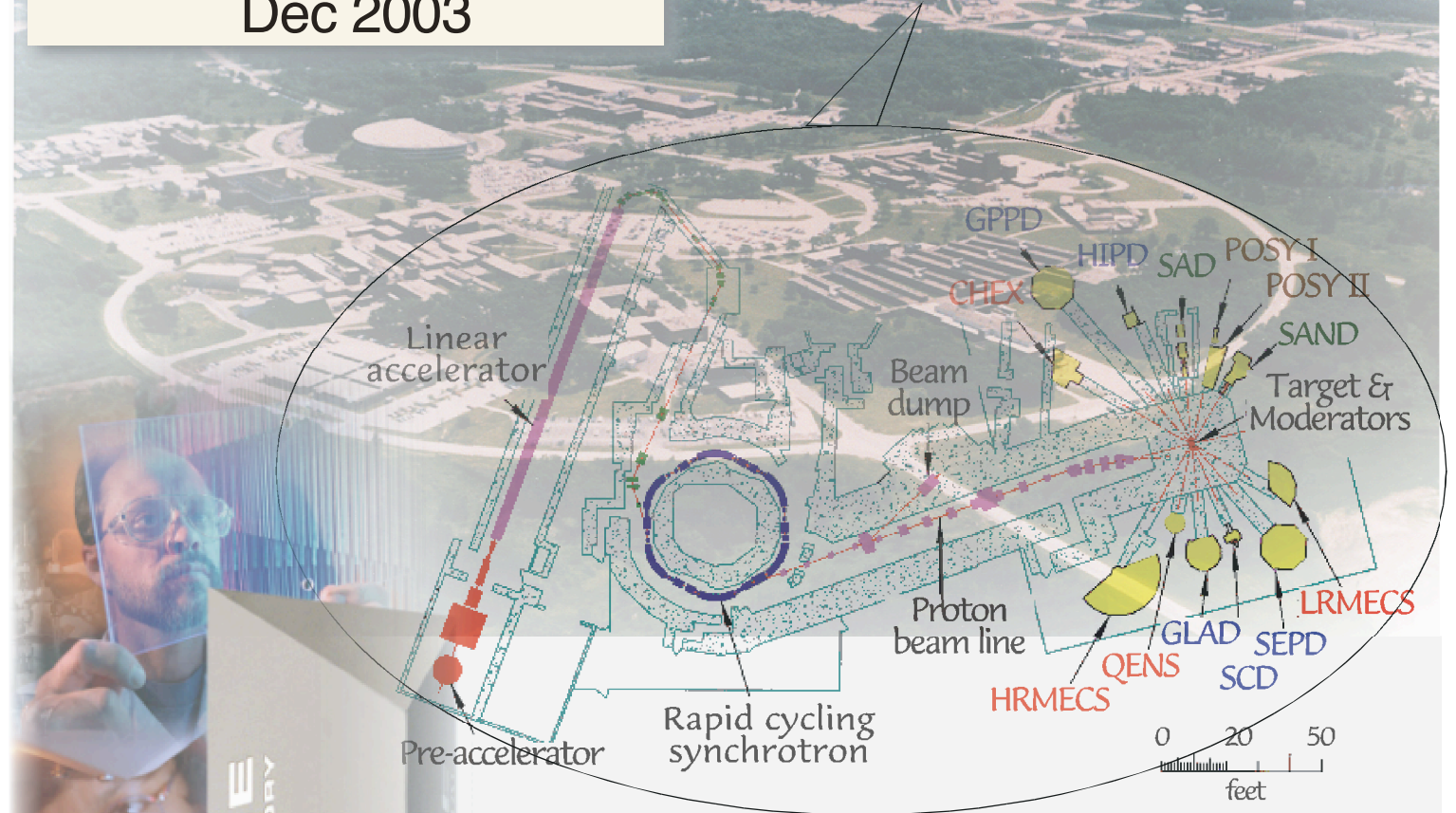

Content

QENS Contacts	3
A Brief History	4
Instrument	
Layout	5
Analyzer-Detector Configuration	6
(Q,E)-Coverage & Energy Resolution	7
Sample Environment	8
Data Acquisition System	10
Data Reduction & Analysis	
QENS Data Reduction & Analysis	11
Check Inelastic Data	12
Calibration and Normalization of Inelastic Data	13
Convert, Integrate, and Overplot Inelastic Datasets	14
Check Diffraction Data	15
View & Overplot Diffraction Data	16
ASCII Files of Diffraction & Inelastic Data	17

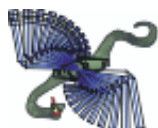

Welcome to the QuasiElastic Neutron Spectrometer (QENS)

This user manual is meant to assist you - the QENS valued customers - to conduct experiments on QENS. We hope to provide the accurate and current information. You are the best person to decide if the manual is useful or needs improvements. Please contact us.

QENS Contacts


Scientific Assistant: Rob Connatser, Rm A235 Building 360
1-630-252-5762 (Office, dial the last 5 digits within ANL)
_____ (cell)
_____ (Home, dial 7 first for outside calls)
Email: rconnatser@anl.gov

Instrument Scientist: Chun Loong, Rm C225 Building 360
1-630-252-5596 (Office)
_____ (Home)
Email: ckloong@anl.gov

Former Visiting Scientist: Jean-Marc Zanotti
Active collaborator on QENS projects Email: jmzanotti@cea.fr

QENS Passwords

Please ask us

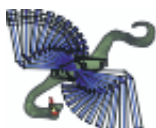

A Brief History of QENS

QENS is the result of a combined effort of many contributors. At the outset it has been a user instrument. The experience of interacting with numerous users for almost two decades has provided valuable guidance for the in-house staff to enhance QENS scientific capabilities. Some key milestones of improvement are:


- 1986 QENS was commissioned following the retirement of an earlier Crystal Analyzer Spectrometer (CAS). It featured three detector arms, housing a total of 35 ^3He detectors. Only a small portion of the scattering angles was covered at each setting. Collection of data over the full Q-range required several measurements each conducted at a different scattering angle setting by rotating the entire detector-arm unit.
- 2001 QENS was completely rebuilt and commenced user operation. It includes 22 crystal-analyzer-detector arms outfitted with cooled beryllium-filters in a new scattering chamber. An average energy resolution (FWHM) varying from 80 meV (elastic) to 4-5% of the energy transfer (inelastic, up to ~200 meV neutron energy loss) was realized over a wide range of wavevectors from 0.3 to 2.5 \AA^{-1} . Two additional detector clusters enabled concurrent measurements of diffraction patterns from 0.1 to 30 \AA^{-1} . Because of the substantial reconstruction resulting in significant enhancement, we often refer the present instrument as the “New QENS”.
- 2003 A supermirror, funnel guide was installed in the incident flight path upstream of the sample, resulted in a neutron-flux gain factor of 2.5 (for 5 \AA -neutrons at sample position).

References:


1. K. F. Bradley, S.-H. Chen, T. O. Brun, R. Kleb, W. A. Loomis, and J. M. Newsam, *Nucl. Instr. Methods*, **A270**, 78 (1988).
2. R. K. Crawford, F. Trouw, and H. A. Belch, IPNS Notes No. 79 (1997), unpublished.
3. R. W. Connatser, Jr., H. Belch, L. Jirik, D. J. Leach, F. R. Trouw, J.-M. Zanotti, Y. Ren, R. K. Crawford, J. M. Carpenter, D. L. Price, C.-K. Loong, J. P. Hodges, and K. W. Herwig, in *16th Meeting of the International Collaboration on Advanced Neutron Sources*, edited by G. Mank and H. Conrad (Forschungszentrum Jülich GmbH, Jülich, Germany, Düsseldorf-Neuss, Germany, 2003), Vol. I, p. 279.
4. QENS web page <http://www.pns.anl.gov/instruments/qens/>



QENS Layout

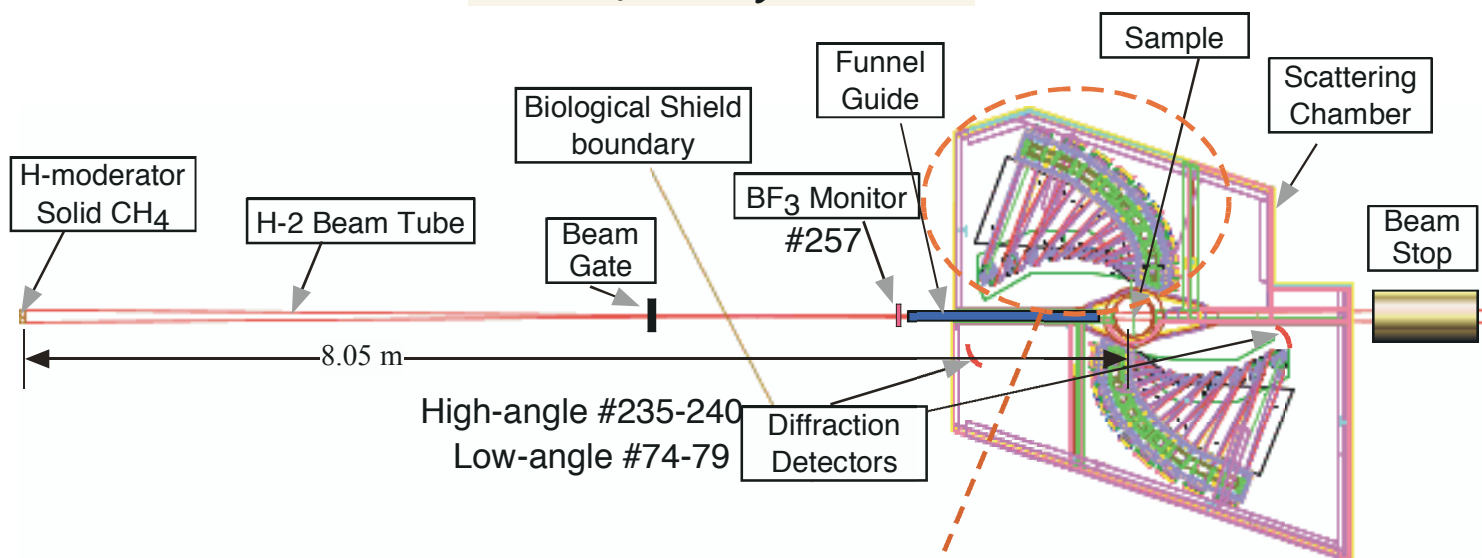

Top view of the analyzer, Be filter, and detector configuration for the high-angle side

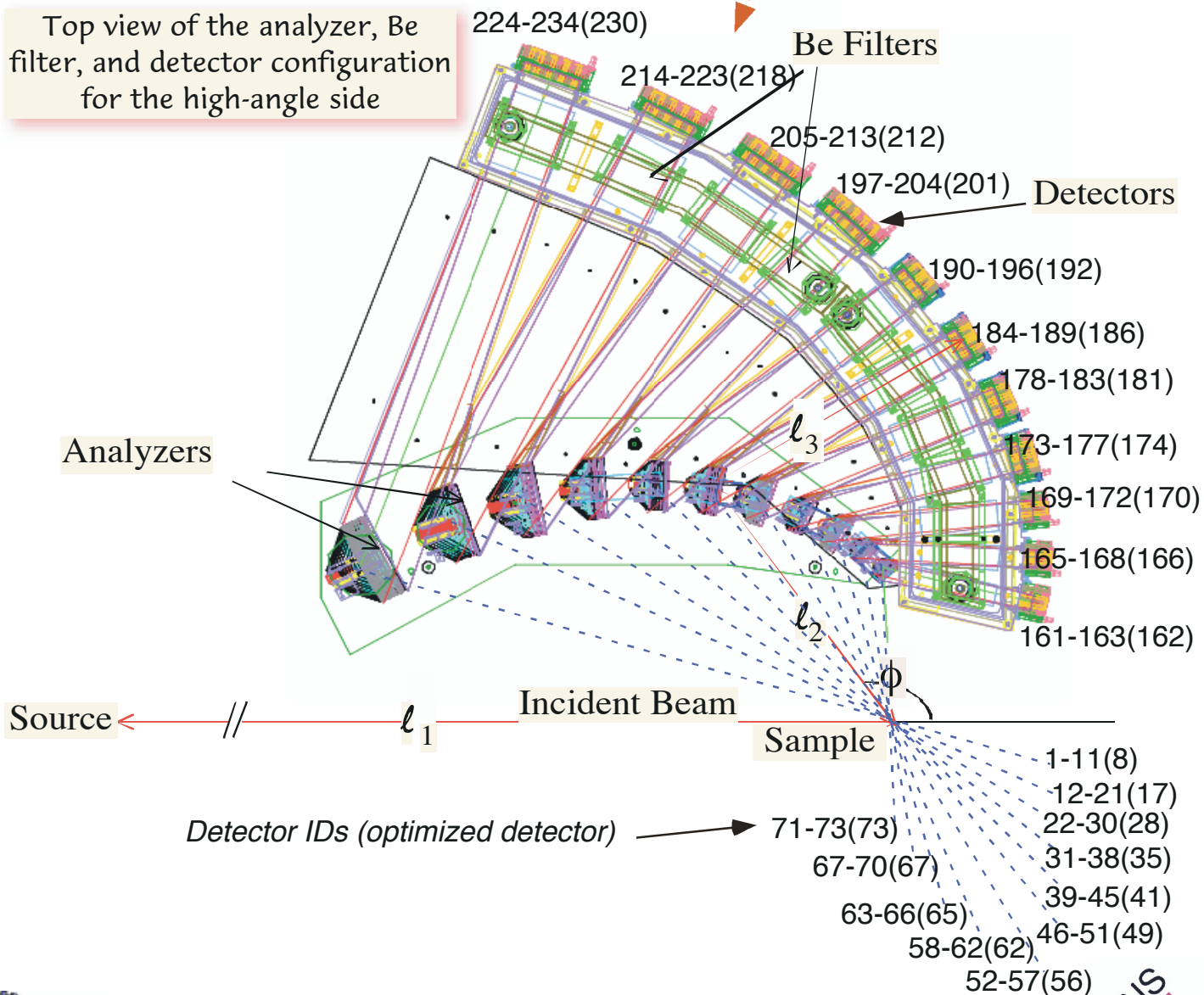

QENS Analyzer-Detector Configuration

Table 1: Parameters for the analyzer-detector arms of QENS



Arm ^{\$}	$\phi(^{\circ})$	$l_2 = l_3$ (m)	$\theta_B (^{\circ})$	E_f (meV)	Analyzer area (cm ²)	No. ³ He counters*
1	93.5	0.61	53.7	2.83	74	3
2	100.5	0.70	52.1	2.95	99	4
3	107.5	0.81	51.0	3.04	167	4
4	114.5	0.92	50.2	3.11	185	5
5	121.5	1.05	49.6	3.17	271	6
6	128.5	1.19	49.0	3.23	296	6
7	135.5	1.35	48.5	3.28	345	7
8	142.5	1.53	48.1	3.32	469	8
9	149.5	1.73	47.7	3.37	629	9
10	156.5	1.94	47.3	3.41	740	10
11	163.5	2.19	47.0	3.44	907	11

*Cylindrical tube: ID = 6.4mm, height = 12.7cm, He pressure = 10 atm.


^{\$} Same parameters for corresponding arms at scattering angles $180^{\circ}-\phi$.

Energy transfer E as a function of flight time at detector t :
$$E = E_i - E_f = E_f \left(\frac{l_1}{l_2 + l_3} \right)^2 \left(\frac{t}{l_2 + l_3} \sqrt{\frac{2E_f}{m_n}} - 1 \right)^{-2} - E_f$$



Note that the TOF at a fixed energy, say the elastic channel, differ for different detector arms. Therefore, it is imperative to measure a vanadium standard for each experiment. The vanadium data will be used to calibrate the energy scale for all the detectors and to normalize the spectra against the incident flux distribution. Do make sure the diffraction detectors are working properly.

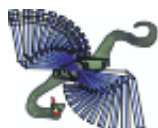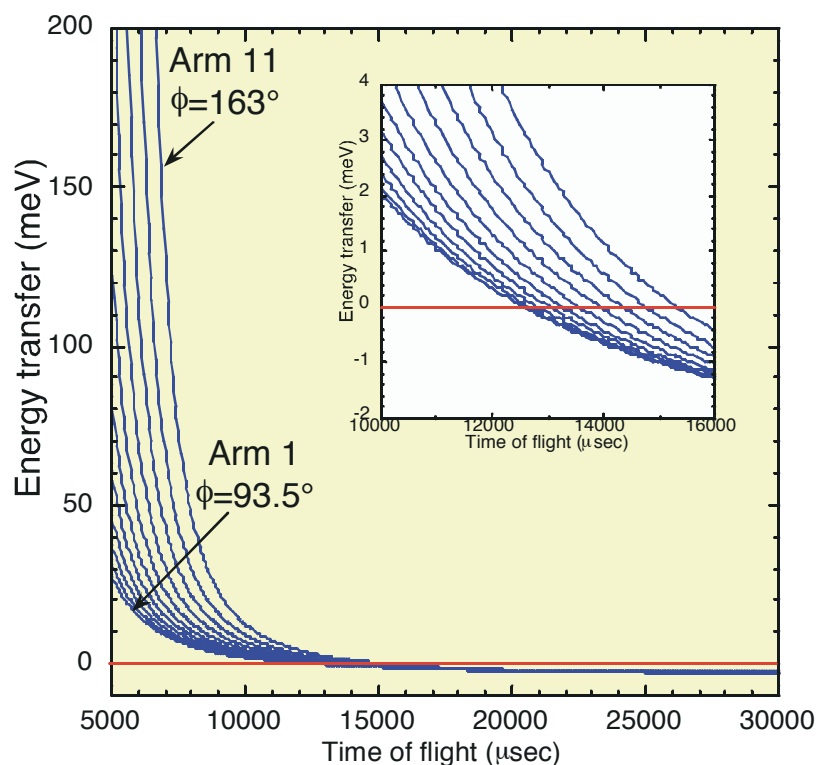

(Q,E)-Coverage

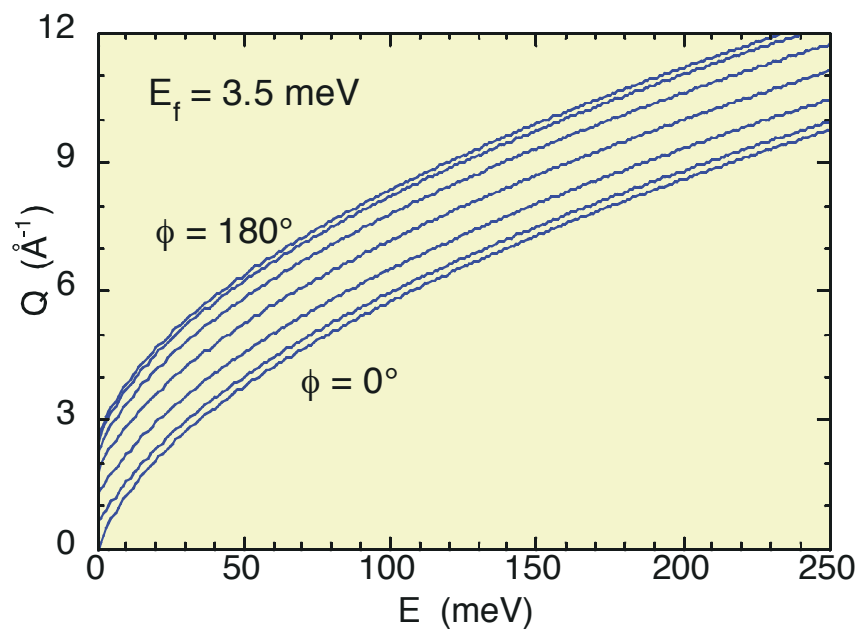

Energy Resolution

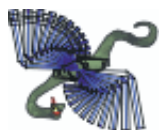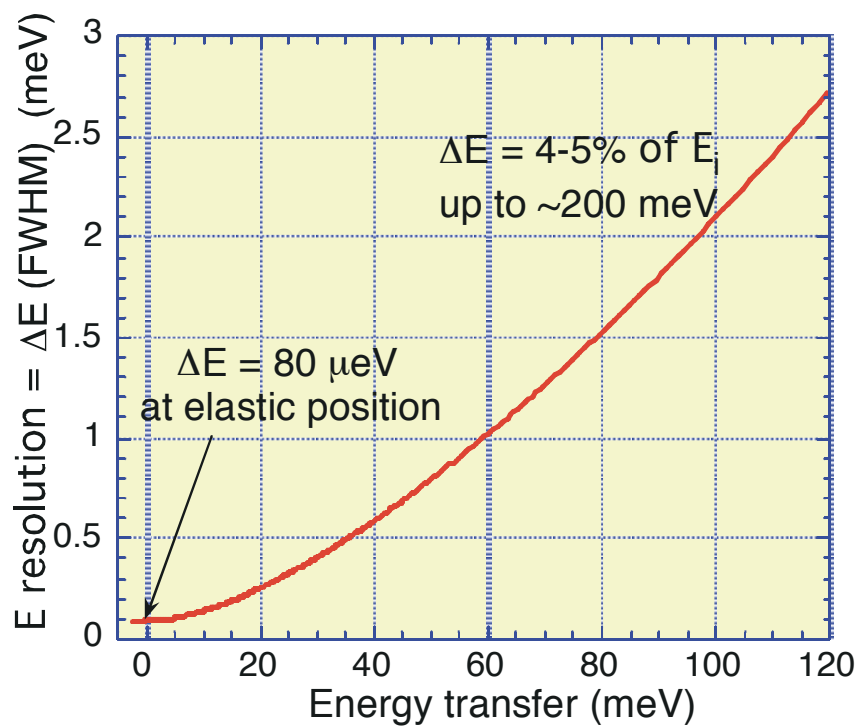

Sample environment

The beam size defined by the funnel guide at the sample position is 3 inches tall by 0.5 inches wide (76mm x 12.5mm). The beam center is 12 inches below the vacuum seal O-ring surface of the sample well. Currently, we are using the sample cans made for the prior beam height of 4 inches.

Ancillary equipments for computer-controlled sample environments:

- Room Temperature Mount
- Instrument Displex: 10-300K
- Hot Stage Displex II: 20K-800K
- Orange Cryostat (50mm diameter bore): 1.5K - 300K
- Orange Cryostat (100mm diameter bore): 1.5K - 300K
- Howe Furnace: RT- 1100°C in vacuum
- Gas handling system for condensing gas molecules in samples

Sample geometry: cylindrical, annular, or flat-plate type. Sample cells are usually made of aluminum alloy and sealed at one end with an indium ring.




Cylindrical cells of id's of 0.125, 0.25 and 0.6 inches are available. Shown here are some of the cylindrical cells for the Displexes and cryostats. *In-situ* gas loading is possible for cylindrical cells.






The annular cell is formed by fitting an inner cylinder into a larger one. Cells of total thickness in the beam of 0.1, 1.0, 1.5 and 2.0 mm are available for use of a Displex. Flat-plate cells of thickness of 0.2, 1.0 and 2.0 mm are available for the Displexes and cryostats.

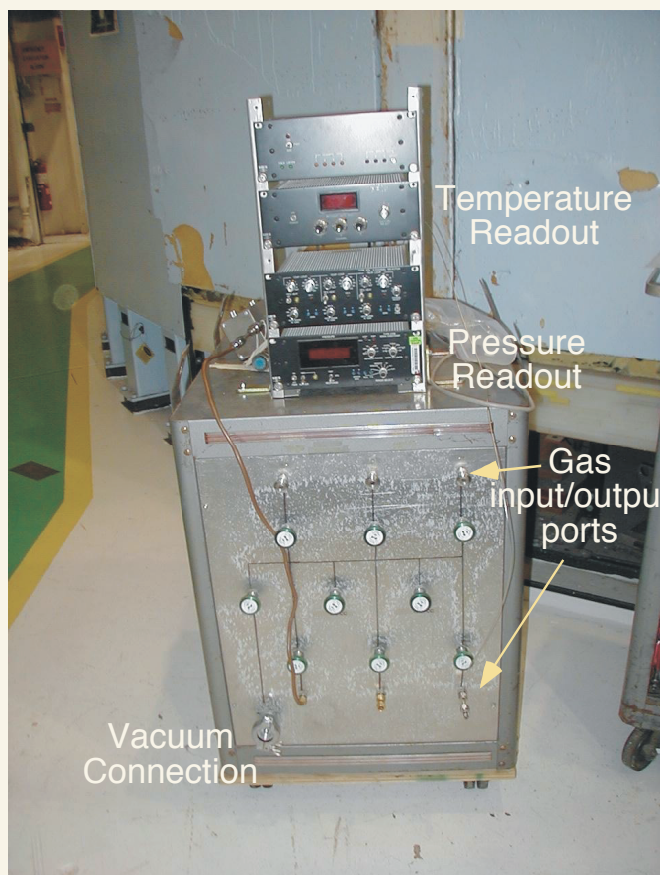

Typically, gas is delivered via a capillary line from the gas handling system through the Displex mounting plate to the sample cell. Prior to gas loading the system is to be leak checked and pumped out thoroughly. The gas from a known volume at a specific pressure and temperature is expanded into the delivery system and absorbed into the sample while cooling.

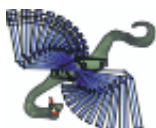


Data Acquisition System

At present QENS is operated using a VAX station as the front-end computer under the IPNS first-generation data acquisition system (DAS).

Communication with the VAX is usually conducted via a terminal-emulation software (e.g, VersaTerm or Exceed) or by Telnet through the local network. Run setup and sample environment control, etc. are entered by line commands. Live run and data archive are viewed in a graphic window one detector at a time. For example, `dis` displays the spectrum of a detector and `cv` clears the display. Do not leave the displays on unattended because over a long period of time it may tie up the computer system and disable all communications. Always clear the video screen display by entering `cv` before leaving the terminal.

QENS is scheduled for an upgrade to new electronics and DAS in the near future. The new DAS - *Integrated Spectral Analysis Workbench (ISAW)* - is a Java-based, platform-independent, Internet accessible, fully portable and highly interactive software package which has been successfully in operation on some of the IPNS modernized instruments. The new hardware will enhance the performance and capabilities of the detectors, sample environment control, and execution of runs.

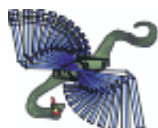

QENS Data Reduction & Analysis


The QENS data treatment software is based on IDL which is available via a multi-user license for the local IPNS networked PCs. Look for the QENS User PC in the User Room which should have the necessary data mapping and preferences set up already for analysis of QENS data.

Since the number of concurrent IDL users is limited, please do not leave IDL idle on the computer. Exit IDL after use.


First, FTP the raw data qensxxxx.run run files from the QENS front-end computer to the PC using WS_FTP95LE (a shortcut on the desktop).

Type a starting string (e.g, qens778* if you want to transfer qens7785.run), then click the Data folder. This will enable a quick look up of the correct runfiles

Data are store in Y:\Zanotti\Qensdata

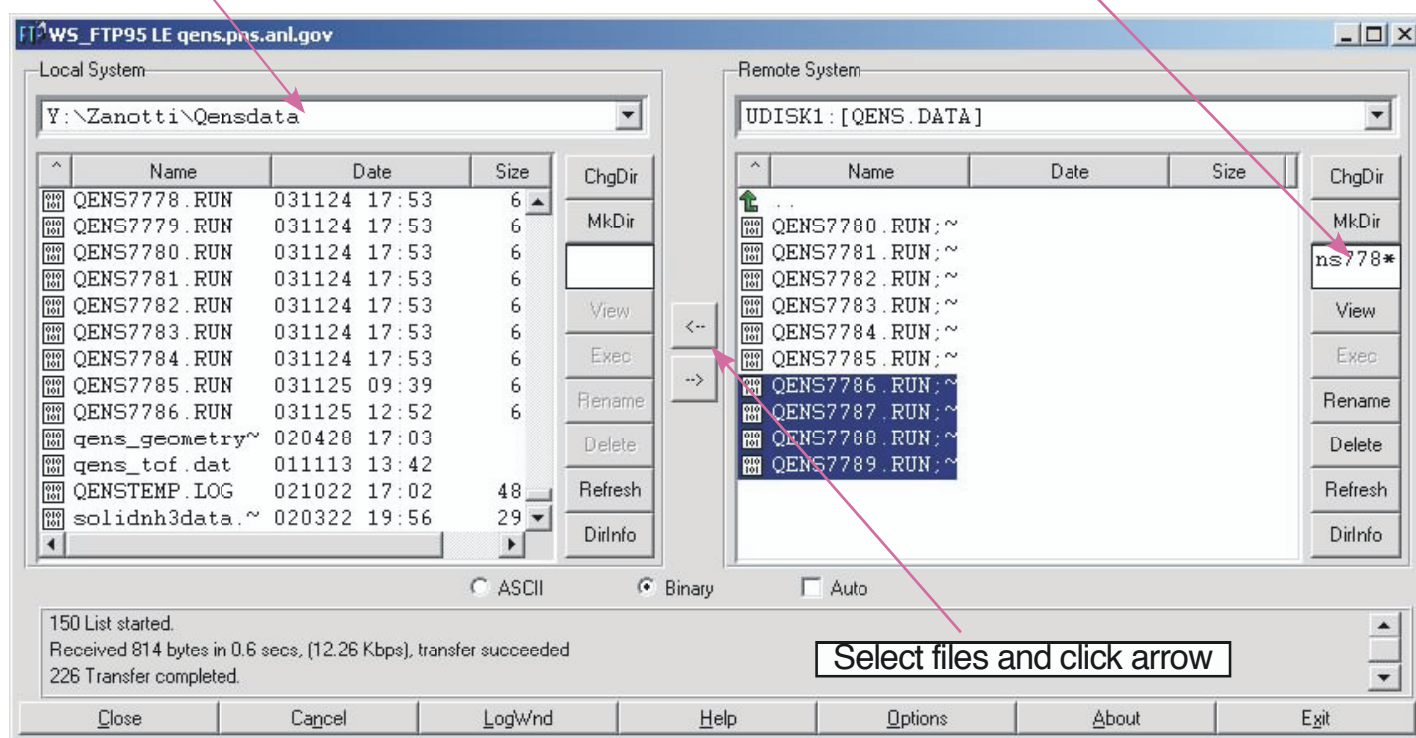


To start IDL, double click the IDL logo, type qensh in the command window.

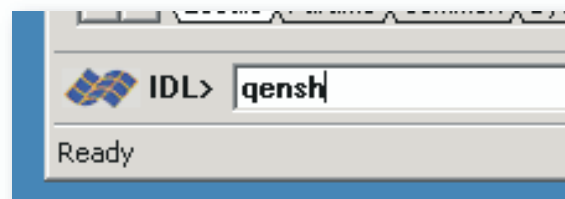

Select **Inelastic>Read Raw data** from the Qensh General Menu. Enter the run number(s) of sample run, select normalization option and an index of the structure (work space) for the data storage in memory.

Check Inelastic Data

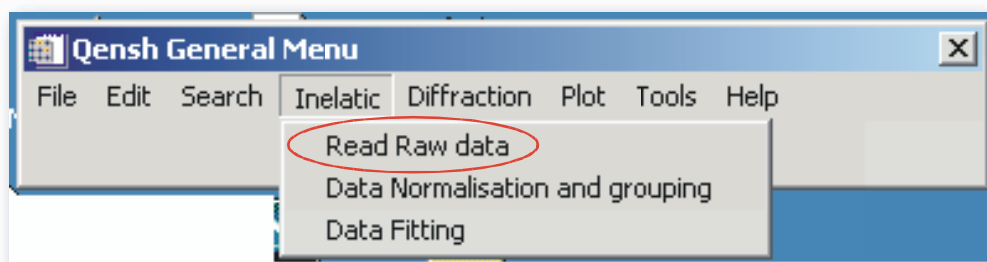

run numbers separated by commas


empty run data will be subtracted from the sample run

options:
time
monitor
pulses
protons on target

run & spectrum information of


The 'Read Raw Inelastic Data' dialog box is shown with the following components and annotations:

- Sample(s) Run #**: Input field containing '5627'.
- Empty Can Run #**: Input field.
- Normalization**: Dropdown menu set to 'Monitor'.
- Read it in structure**: Button.
- work space index**: Input field containing '2'.
- Current IDL index #**: Input field containing '84'.
- IDL index, not detector ID**: Annotation pointing to the 'IDL index #' label on the left plot's y-axis.
- Left Plot**: A scatter plot of 'IDL index #' (y-axis, 0 to 150) vs 'tof (microsec)' (x-axis, 0 to 2.0x10⁴). It shows a dense cluster of points with a red horizontal line at IDL index ~84.
- adjust z-axis max & # contours**: Annotation pointing to the '# of levels = 10' and 'Zmax = 200' controls at the bottom left.
- double click to select a detector for display on the right**: Annotation pointing to the '<<' and '>>' buttons below the left plot.
- Right Plot**: A line plot of 'Intensity (A.U.)' (y-axis, 0 to 400) vs 'tof (microsec)' (x-axis, 1.00x10⁴ to 1.60x10⁴). It shows a sharp peak at approximately 1.18x10⁴ microseconds. Text above the plot reads: 'current_WS = 2', 'Q0=2.05 Å⁻¹', 'Norm.=Monitor', 'Detector ID# 173', and '2θ=114.°'. Below the plot, it says 'YBC 0.2 N2+moist 200C after heating to 700C'.
- cursor position**: Annotation pointing to the 'X = 12780.0' and 'Y = -312.585' readouts.
- x- & y-axis limits**: Annotation pointing to the 'Xmin = 10000', 'Xmax = 15000', 'Ymin = 0.', and 'Ymax = 400' controls.
- error bar option**: Annotation pointing to the 'Error Bars' checkbox, which is checked.
- scroll through detector IDs**: Annotation pointing to the '<<' and '>>' buttons below the right plot.

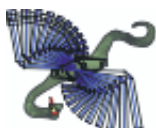

Calibration & Normalization of Inelastic Data

Select **Inelastic>Data Normalization and grouping** from the Qesh General Menu.

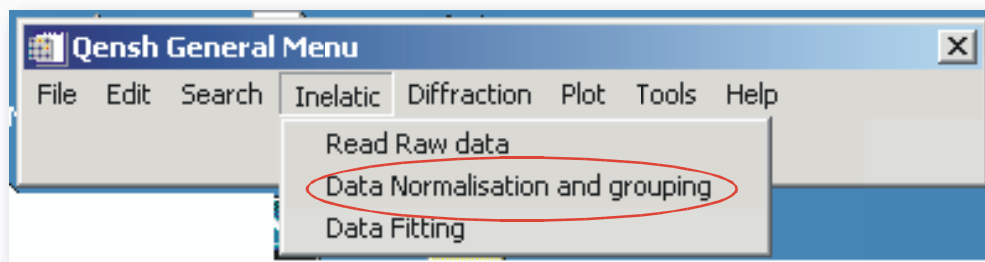

You may use this window many times to generate energy-calibrated, detector-efficiency-normalized datasets in different memory workspace for further comparison or treatment later.

run numbers separated by commas; empty run(s), if given, will be subtracted from data

vanadium run(s), if given, will be used to normalize data according to detector efficiencies

The screenshot shows the 'Data Normalization and grouping' window. On the left, there are input fields for 'Sample run(s)' (6563), 'Cell for Sample', 'Vanadium', 'Cell for Vana.', 'Buffer', and 'Cell for buffer'. Below these are buttons for 'Read it in structure' and 'Group it in Structure'. A 'Bad detectors list' is shown with the text '11,21,31,63,68,69,70,165,167,168,169,170,171,172,214,10,'. On the right, there is a plot titled 'e Furnace, SS can, empty with normal a' showing 'Intensity (A.U.)' vs. a scale from -1.0 to 2.0. The plot shows a sharp peak at 0.0. To the right of the plot are input fields for 'Xmin = -1.4000', 'Xmax = 2.00000', 'Ymin = 0.', and 'Ymax = 0.5e5'. At the bottom, there are buttons for '<<', '10', '>>', and 'Plot it'.

for subtraction of buffer component from sample

x- & y-axis limits setting

Generate spectra from individual detectors in memory workspace

Generate spectra grouped according to detector arms in memory workspace


Scroll spectra through detectors or groups

May edit bad-detector list to eliminate them from data




Convert, Integrate, and Overplot Inelastic Datasets

Select **Plot>Inelastic** already grouped... from the Qesh General Menu.

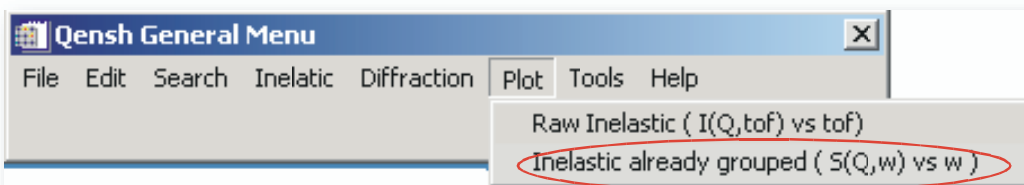


This GUI displays the previously generated datasets from memory workspace as the full structure factor $S(\phi, E)$ (ϕ is the scattering angle) or the vibrational density of states $g(\phi, E)$, or $S(E)$ and $g(E)$ by summing over the detector groups.

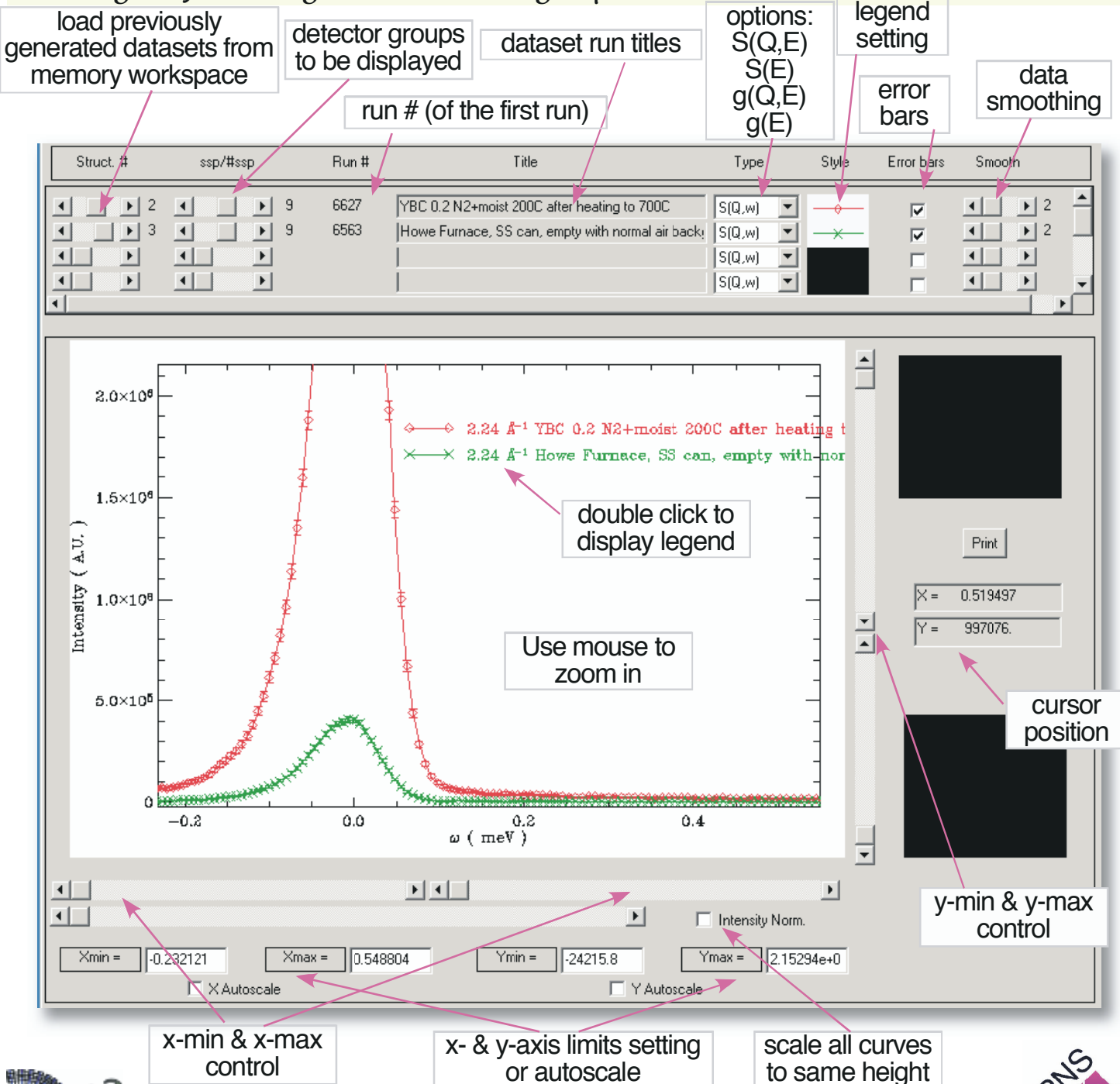


Check Diffraction Data

Select **diffraction** > **Read Data** from the Qensh General Menu. Enter the run number and an index of the structure (work space) for the diffraction data storage in memory. Click the "Read in Structure #" box.

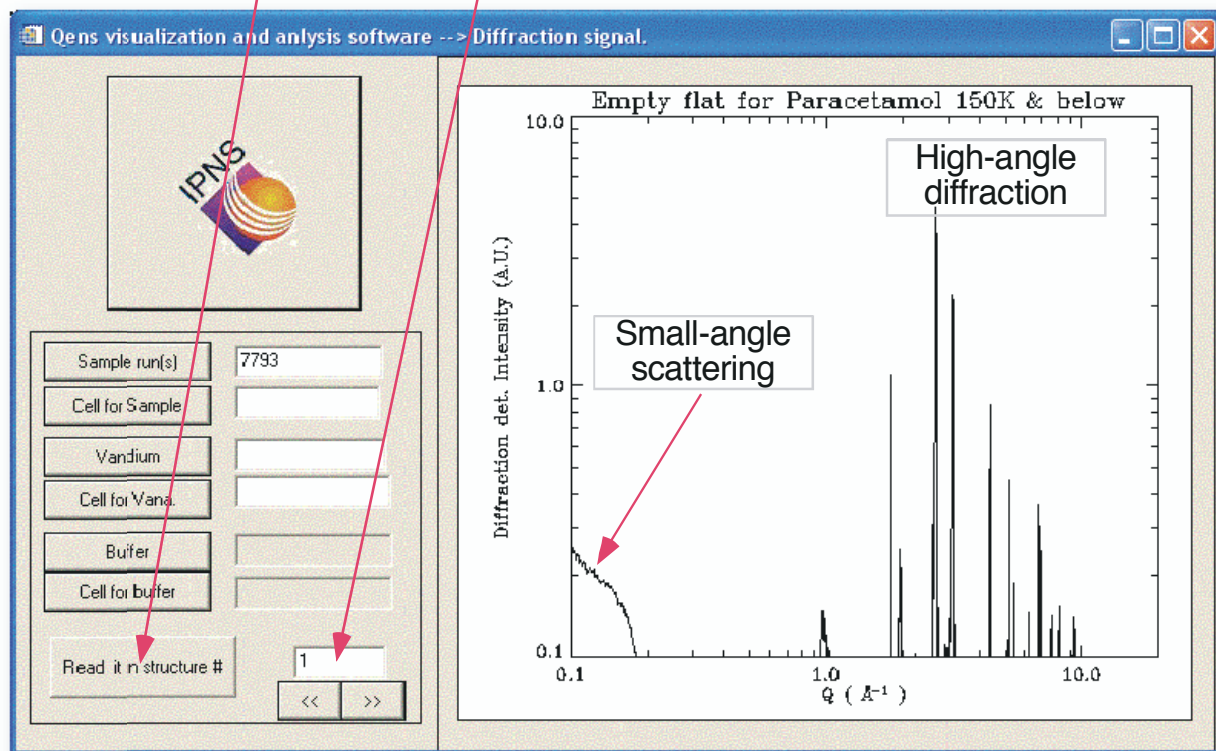


The spike seen here is an artifact. These plots in log-log scale is to show an overview of the diffraction data. For detailed plot, choose **Plot > Inelastic... S(Q,w) vs q** as shown below

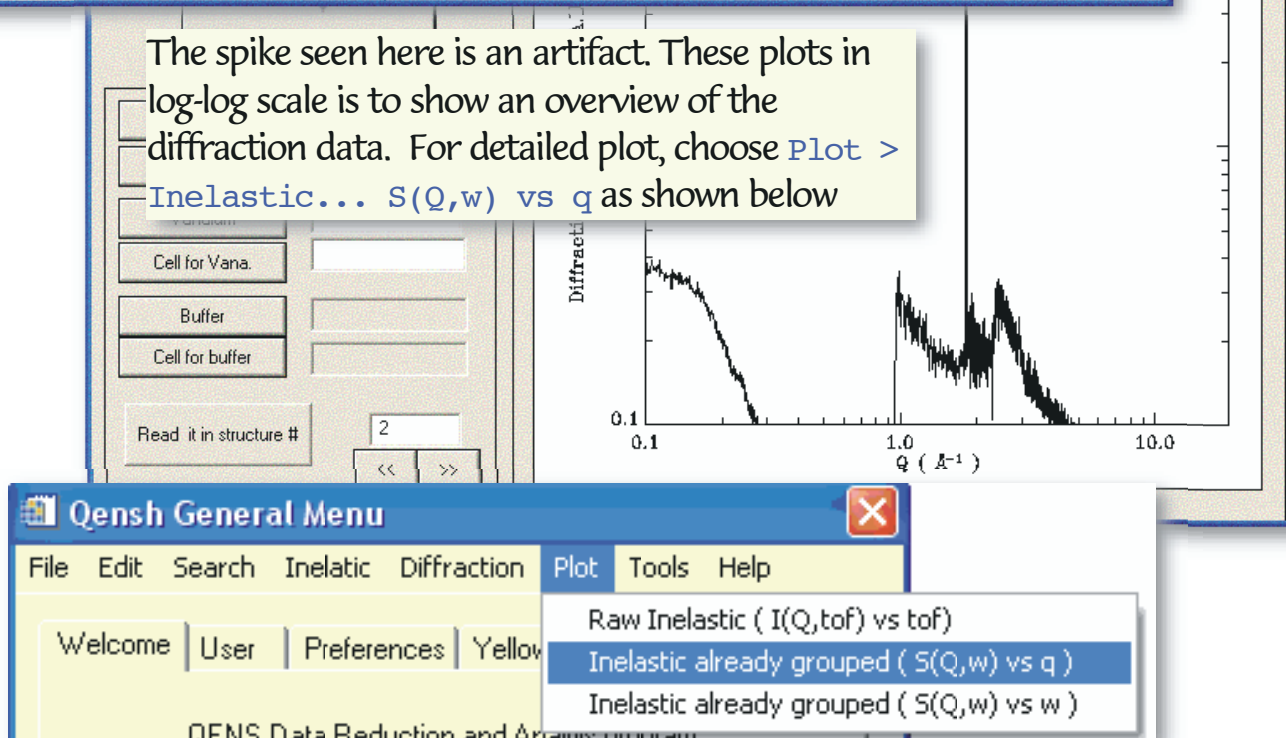

View and Overplot Diffraction Data


Previously generated
diffraction data sets from
memory workspace

detector index

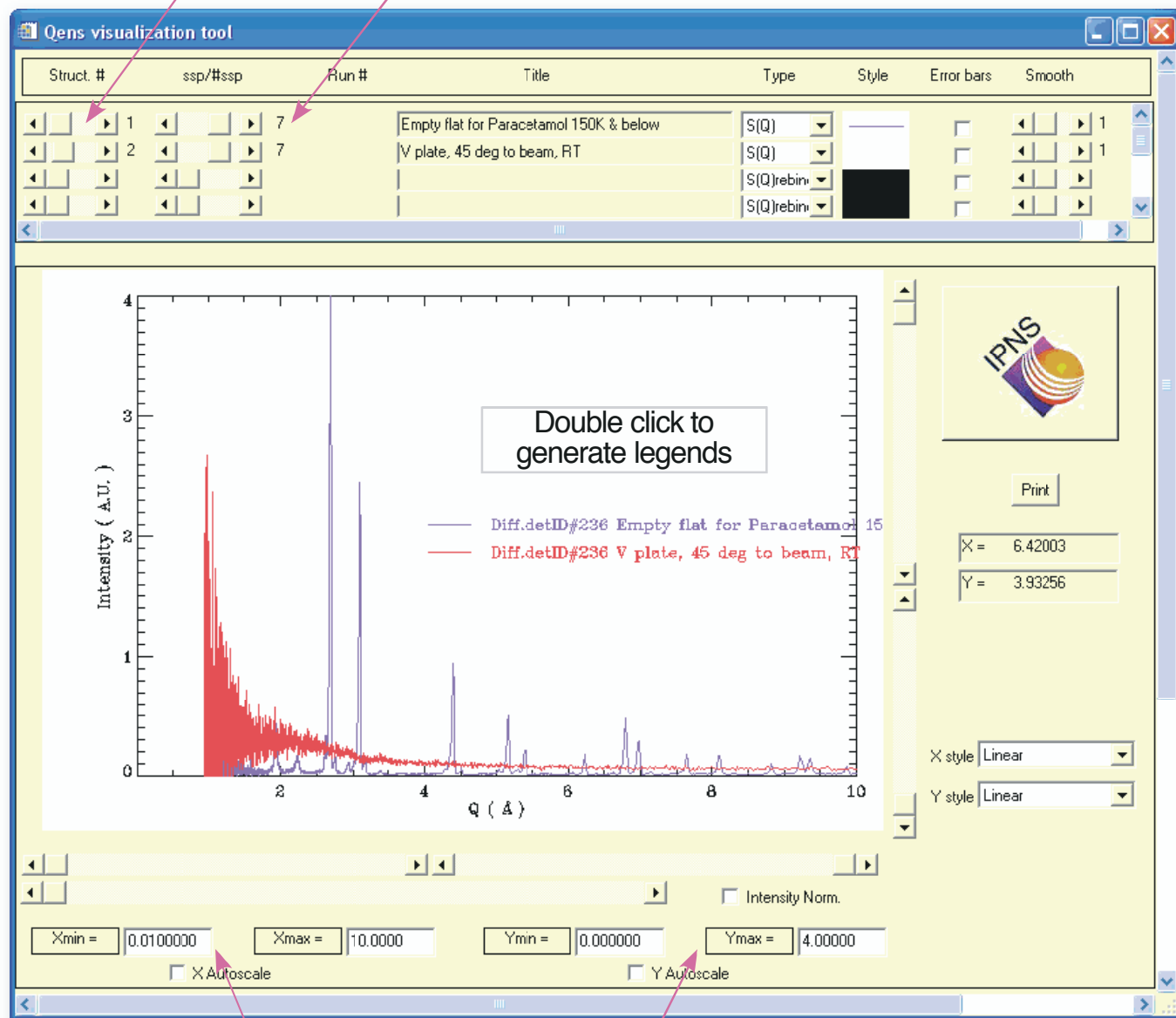

x-min & x-max control

x- & y-axis limits setting
or autoscale

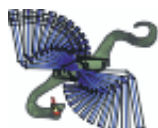

ASCII Files of Diffraction & Inelastic Data

ASCII files of $S(Q)$, $S(\text{detector arm}, E)$, $S_{\text{sum}}(E)$ and $G_{\text{sum}}(E)$ are automatically generated whenever `Plot>...` is executed. The files are stored in `Y:\Zanotti\QENS\datatreatment\`. For example:

Name	Size	Type
Qens_diffraction_7749-n-.dat	891 KB	DAT F
Qens_diffraction_7793-n-.dat	890 KB	DAT F
Qens_u2_7743-n-.dat	1 KB	DAT F
Qens_SQW_7743-n-.dat	821 KB	DAT F
Qens_gw_7743-n-.dat	167 KB	DAT F
Qens_elast_int_7743-n-.dat	2 KB	DAT F
Qens_u2_7743-n7743-.dat	1 KB	DAT F
Qens_SQW_7743-n7743-.dat	821 KB	DAT F
Qens_gw_7743-n7743-.dat	167 KB	DAT F
Qens_elast_int_7743-n7743-.dat	2 KB	DAT F
Qens_eff_7743.dat	14 KB	DAT F
Qens_SQW_7793-n-.dat	821 KB	DAT F
Qens u2 7793-n-.dat	1 KB	DAT F

Data columns: For diffraction:

`q_4345, dspacing_4345, sq_4345, esq_4345, q_sa_4345, d_sa_4345, sq_sa_4345, sqsa_ns5_4345, sqsa_ns10_4345, q_ha_4345, d_ha_4345, sq_ha_4345, sqha_ns5_4345, sqha_ns10_4345, qb_4345, db_4345, sqb_4345`

correspond to: Q (\AA^{-1}), d-spacing (\AA), $S(Q)$, $S(Q)$ -error, ...

and sa/ha designates the sum over the small/high-angle detectors

nsx, ... are data smooth over x points

qb, ... are rebinned data of all detectors according to the common Q values.

For inelastic data:

`w (meV), d_0_0.36_A-1, e_0_0.36_A-1, d_1_0.36_A-1, e_1_0.36_A-1 ...`

correspond to E (meV), $S(E)$ for det group 0 with elastic Q of 0.36 \AA^{-1} , $S(E)$ -err..., etc.

